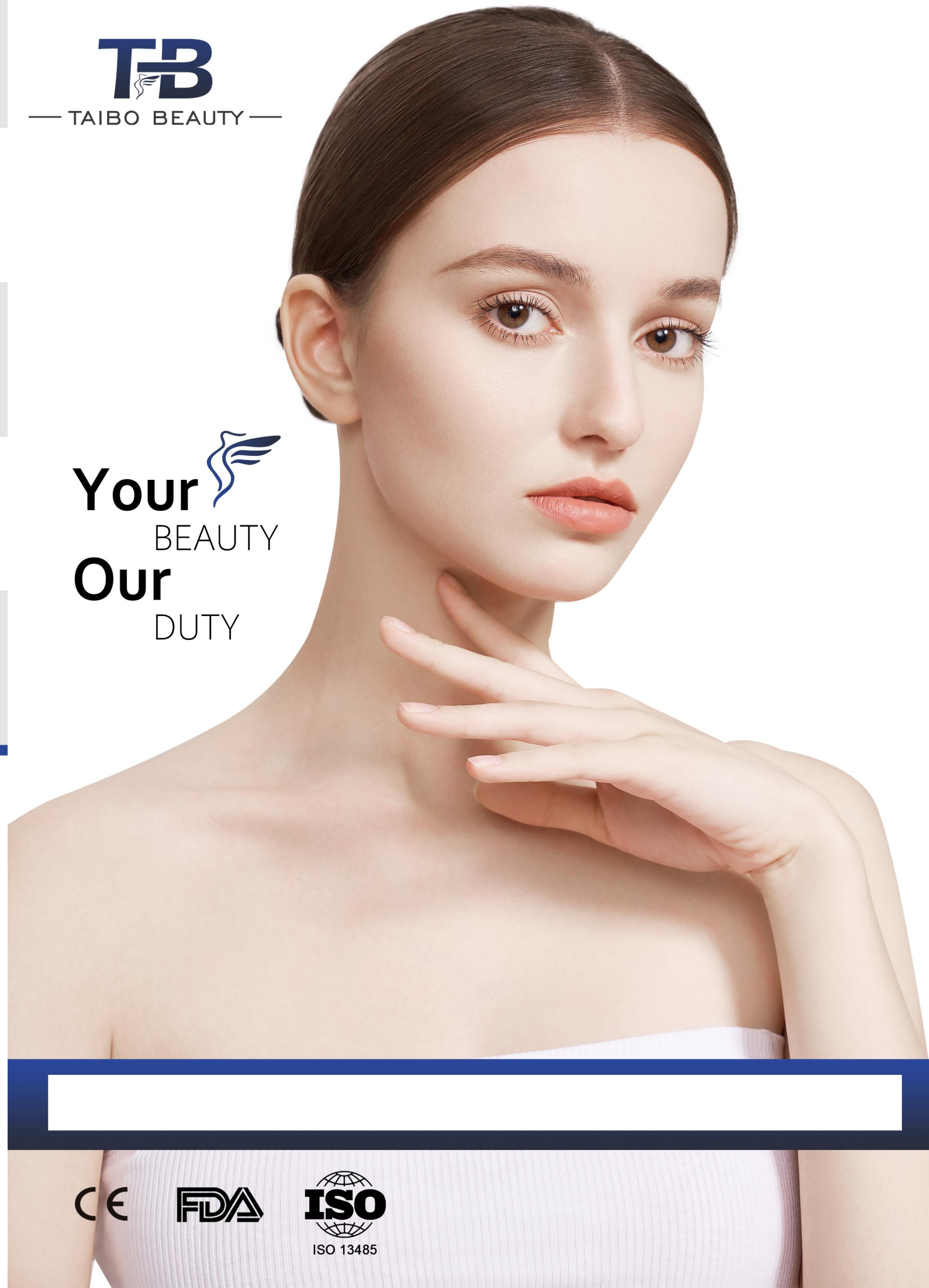




TB
— TAIBO BEAUTY —

Your 
BEAUTY
Our
DUTY



CE

FDA

ISO
ISO 13485

DepiLazer[®]

808/808+755+1064+940NM
DIODE LASER MACHINE

Imported laser (US Coherent) | 15.6-inch large display | Advanced water pump



Superiority

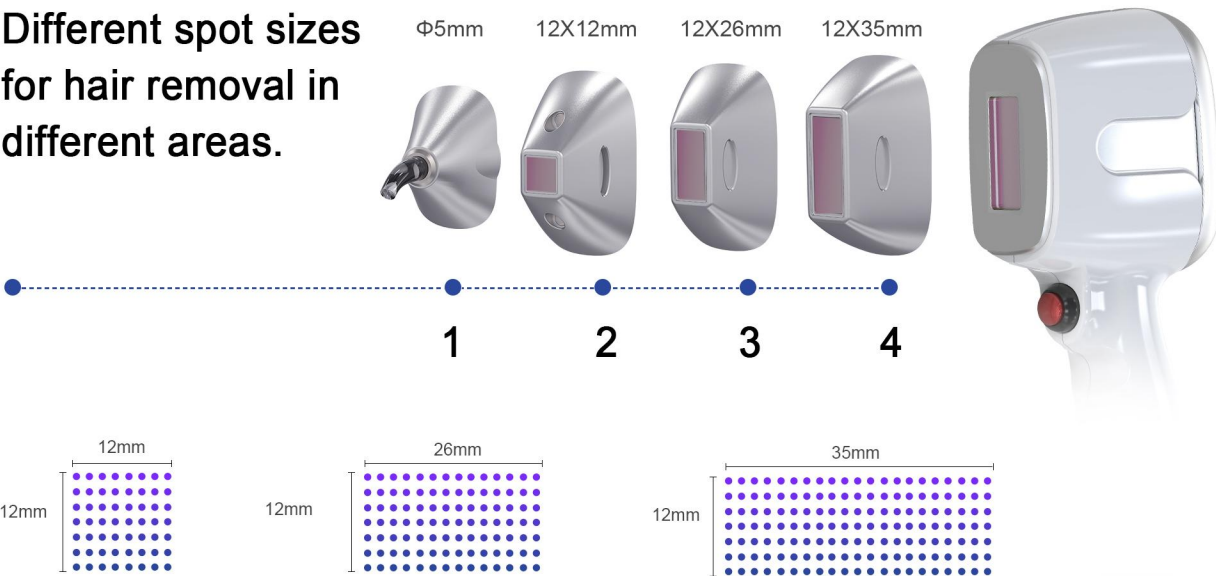
- ◆ Max 50 millions shots
- ◆ 48 hours continuous
- ◆ Collimated Laser
- ◆ Handle Power:2000W

Working principle

The hair is trimmed and a cooling get applied to the area to be treated In a fraction of second pulsed red light heats up the hair The hair becomes detached from its follicle and eventually falls out The hair follicle is left unable to produce a new hair.

Magnetic Head Series

Different spot sizes
for hair removal in
different areas.



Specification

product name	Semiconductor laser hair removal device	Ambient operating temperature	10℃~30℃
Model	TB-808X	Relative working humidity	≤80%
Screen size	15.6 inch	Atmospheric working pressure	700hPa~1060hPa
Treatment heads quantity	1 piece / 2 piece	Storage temperature	10℃~40℃
Output Power	800W/1200W/1600W/2000W(laser power)	Storage relative humidity	≤80%
Maximum input power	5600W	Storage atmospheric pressure	700hPa~1060hPa
Frequency	1-10Hz adjustable	Storage location	0℃~30℃
Wavelength	808nm+940nm+755nm+1064nm (support customization or any combination of wavelength)	Spot size	12X12mm/12X26mm/12X35mm
Cooling System	Air cooling + Water cooling+Electronic cooling +Semiconductor cooling+Compressor cooling	Maximum storage time	Long term
		Packaging material	Air box and foam
		Electrical classification	AC 220V-240V~, 50Hz/ AC 100-120V~, 60Hz

DepiLazer®

808/808+755+1064+940NM
DIODE LASER MACHINE(1600/2000W)

Imported laser (US Coherent) | 15.6-inch large display | Advanced water pump



Superiority

- ◆ Painless hair removal
- ◆ Skin rejuvenation
- ◆ **12x36mm** Spot size
- ◆ **50** millions shots
- ◆ Five cooling modes



4 in 1 PERSONALIZED PLATFORM

808nm,755nm,940nm,1064nm wavelength

From single wavelength to 4 wavelengths, we can customize whatever you want.



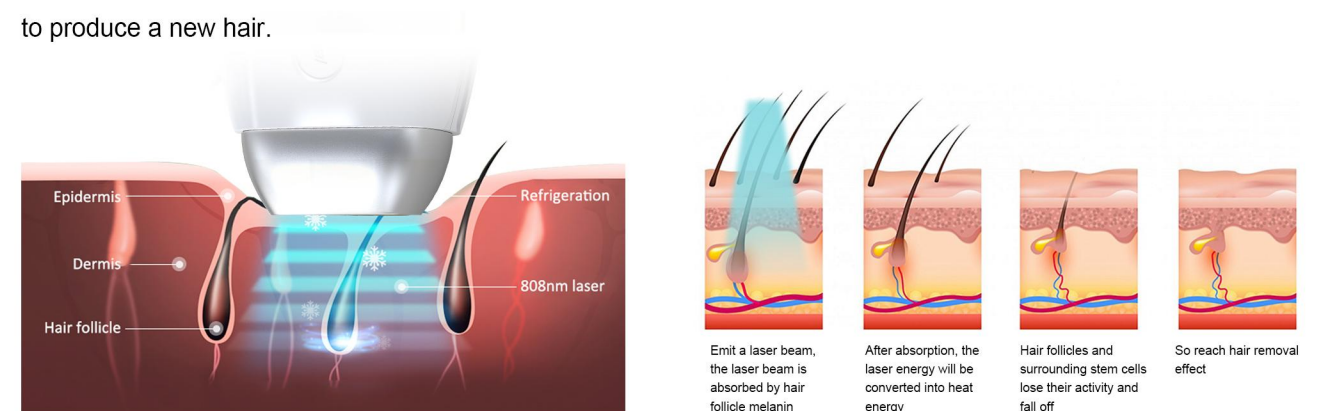
1 minute Cooling **-29°C**

Best **TEC+Compressor Cooling**
system hair remover laser



Working principle

The hair is trimmed and a cooling gel applied to the area to be treated. In a fraction of a second, pulsed red light heats up the hair. The hair becomes detached from its follicle and eventually falls out. The hair follicle is left unable to produce a new hair.



Machine Interface



DepiLazer[®]

808NM/808+755+1064NM DIODE LASER MACHINE



Superiority

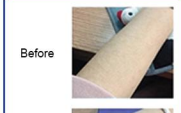

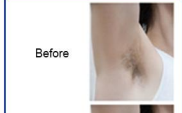

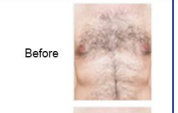
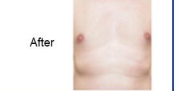
- ◆ Painless hair removal
- ◆ Skin rejuvenation
- ◆ **12x36mm** Spot size
- ◆ **50** millions shots
- ◆ **5** cooling modes



Specification

Laser type	Semiconductor laser	Spot size	12X12 mm
wavelength	808nm	Frequency	1-10HZ
Output Mode	pulse output	Packing size	52X45X66cm
Control Method	Touch control	Net weight	30kg
Screen Size	10.4 Inch	Fuse specification	Ø5×25 10A
Energy Density	2-100J/cm2	Light spot surface temperature	-5—5℃
Pulse Width Range	10-400ms	Power	AC220V±10% 10A 50HZ , 110v±10% 10A 60HZ
Input power	1800w	Cooling System	water cooling+ electronic cooling + semiconductor cooling
output power	300W		

Treatment parameters

<p>Information</p> <ul style="list-style-type: none">• Gender: Female• Age:27• Country: China• Skin Type: Yellow skin <p>< 808nm hair removal ></p> <ul style="list-style-type: none">• Energy(J/cm²): 80• ePulse Width(ms): 150ms• Frequency: 1-5	<p>Before</p>  <p>After</p> 	<p>Information</p> <ul style="list-style-type: none">• Gender: Female• Age:31• Country: Sweden• Skin Type: Very white skin <p>< 808nm hair removal ></p> <ul style="list-style-type: none">• Energy(J/cm²): 95• ePulse Width(ms): 300ms• Frequency: 1-5	<p>Before</p>  <p>After</p> 	<p>Information</p> <ul style="list-style-type: none">• Gender: Male• Age:25• Country: Canada• Skin Type: White skin <p>< 808nm hair removal ></p> <ul style="list-style-type: none">• Energy(J/cm²): 90• ePulse Width(ms): 200ms• Frequency: 1-6	<p>Before</p>  <p>After</p> 
---	--	--	--	---	--

ALEXANDRITE LASER HAIR REMOVAL
1064NM/755+1064NM

8-22mm to adjust the spot size | Optical fiber imported from the Germany 1.5mm



Superiority

- ◆ Skin Tightening
- ◆ Hair Removal
- ◆ Crystal Skin Cleansing
- ◆ Acne
- ◆ Onychomycosis
- ◆ Vascular skin lesions

Working principle



The gentle hair removal laser targets the hair follicle and thermally destroys it with intense laser light energy, providing permanent hair removal over the course of several treatments.

Magnetic Head Series



Germany imported optical fiber

The transmission of energy through optical fiber is more stable



Infrared Aiming Beam

Make the treatment more accurate



755nm+1064nm

Wavelength



8-22mm

Adjustable Spot Size



0.25-100MS

Pulse Width



Best cooling system

Water&TEC&FAN



Specification



Screen size	15.6 inch	Fiber	Germany Imported 1.5mm
Wavelength	1064nm/755nm	Skin coolingsy system	semiconductor cooling
Pulse width	0.25ms-100ms	Indicator light	Adopts 532nm greenlight
Maximum Output Energy	1064nm:80J (Max 100J) 755nm:50J (Max 70J)	Repetition Rate	0.5Hz-10Hz
Output power	3000W	Cooling System	Water&TEC&FAN
Spot size	8/12/14mm;16/18/22mm	Power Supply	AC220V 10A 50/60Hz

Non-invasive Fat Dissolving | Muscle Building | Private Anti-Age



Working principle

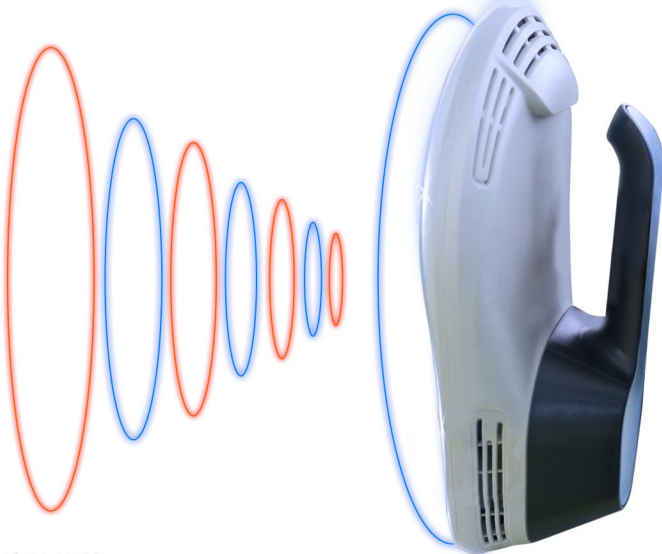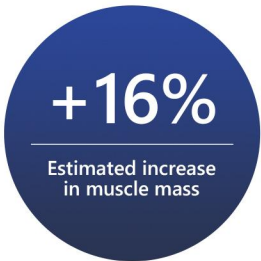EMS uses the latest High-intensity Electromagnetic technology, which has the ability to induce supramaximal muscle contractions. The rapidly changing magnetic field induces electric currents in the tissue where it epolarizes neural membranes and governs motor units in the target muscles, causing concentric contractions.



30 MINUTES LYING DOWN ≈ 5.5 HOURS OF EXERCISE

30 MINUTES ABDOMINAL MUSCLE WORK ≈ 30,000 SIT-UPS

Specification

product type	Pulse electromagnetic wave	Operating Voltage	AC 220V-240V~, 50Hz
Screen	15.6 inch	Total power	6000W
Pulse frequency	1-250HZ	Handle configuration	2 body handles, 2 limb handles
Pulse voltage	500V-1500V	Electromagnetic energy	8 Tesla
Energy ratio	10%-100%	Ambient working temperature	-10℃~30℃
operating hours	1-90min (adjustable)	Relative working humidity	≤80%
Frequency switching	3 kinds of frequency intelligent adjustment switch	Atmospheric pressure range	700hPa~1060hPa
Operating mode	Smart mode	Way of working	Two handles work at the same time/ four handles work at the same time
storage method	Support three modes of data saving		



1-250HZ
1-250HZ ADJUSTABLE



20000 TIMES
20000 TIMES SQUEEZE
MUSCLE IN 30 MIN



4 MODES
4 TREATMENT MODES





EFFICIENT
SPEED UP FAT BURNING

Non-invasive Fat Dissolving | Muscle Building | Private Anti-Age



Application

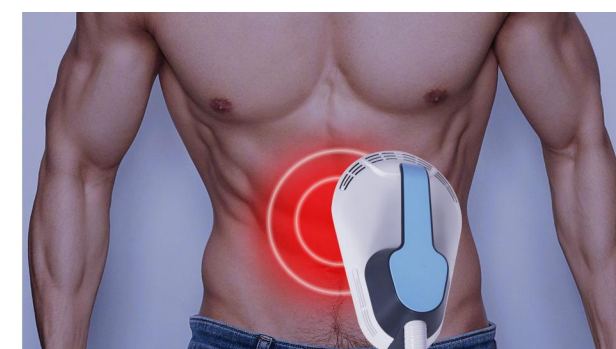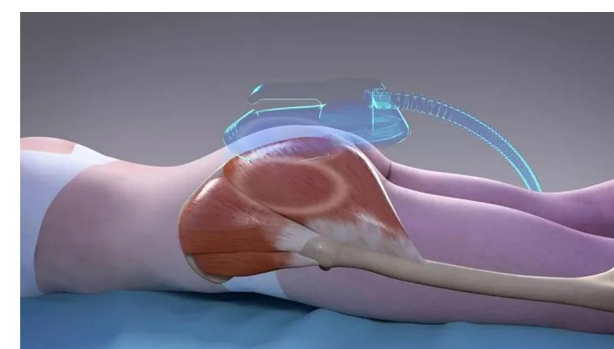- ◆ Hip lifts, waistcoat line
- ◆ Men who need to gain muscle change
- ◆ People who need to lose weight
- ◆ Improving the leakage caused by the relaxation of the basin muscle
- ◆ Postpartum mother(separation of rectus abdominis) - improve the shape of the abdominal muscles and shape a flat abdomen.



21% ↓

EFFECTS ON FAT

Several recent studies using CT MRI and ultrasound evaluations have reported approximately 21% reduction in subcutaneous fat layer in patients treated by based device on their abdomen.

16-18% ↑

EFFECTS ON MUSCLE

Recent studies reported that on average 16-18% increase in abdominal muscle thickness was observed in treated patients one to two month after treatments.

“ SELECTIVE PURCHASE

Women's Privilege -Pelvic Floor Treatment

The high-intensity focused electromagnetic field stimulates motor neurons to control the pelvic floor muscles to perform 20,000 Kegel exercises, which is non-invasive and non-invasive to improve pelvic floor muscle relaxation and other symptoms.



Increase muscle



Non invasive



Reduce fat



Patented technology



Non hyperthermia



Noninvasive hip augmentation



4D SILICON 360
CRYOLIPOLYSIS MACHINE

Up to 4 Areas Treated Simultaneously | Thermal shock lipolysis



Application



Special antifreeze film
Avoid frostbite and protect the subcutaneous organs.



One 'press' and one 'install'
Probes are easy to replace plug-and-play plug-in probes, safe and simple



Comprehensive treatment
360-degree refrigeration without dead corners, larger treatment area



Two handles and probes
Double handles and double heads can work at the same time or independently



Heating mode
A 3-minute heating stage can be selectively performed before cooling to accelerate local blood circulation.

7 Different Cartridges Available



Wing Mini
Double Chin,
Axilla,Knees



Flat Mini
Double Chin,
Axilla,Knees



Flat
Abdomen,Arms,
InnerThighs



Wing
Abdomen,Flanks,
Love Handles



Wing Plus
Abdomen,Bra Line,
Flanks,Love Handles,
elow Buttocks



Wing Max
Abdomen,Flanke
ove Handles



Arc-T
Outer Thighs
nner Thighs

Reincent Fat Freezing Machine

360 Cryolipolysis machine uses state-of-the-art freezing technology removes excess subcutaneous fat from a specific location. Designed with 360° Surround Cooling Technology, the multiple cups can fit all curves and contours from as high as your chin all the way down to the knees, to deliver an absolute non-invasive fat reduction way. The ideal temperature from -5 °C to -11 °C which can induce adipocyte apoptosis is cooling energy to achieve non-invasive and powerful lipidlowering.Different from adipocyte necrosis, adipocyte apoptosis is a natural form of cell death. It's to maintain the stability of the internal environment. Cells die in an autonomous and orderly manner, thereby effectively reducing fat cells without causing damage to surrounding tissues.

Treat 4 Areas Simultaneously in Only 40 Minutes

- ◆ 360 Cryolipolysis machine will come with 4 handles can work at the same time.
- ◆ Each handle is individually controlled and can be adjusted with different parameters.
- ◆ Available to treat multiple areas or multiple people simultaneously.

Specification

Product name	4D silicon 360 cryolipolysis machine	Vacuum suction	1-5 gears (10-90Kpa)
Display screen	10 inch large LCD	Setting time	1-99min (default 60min)
Cooling temperature	1-5 gears (cooling temperature 0 to -16℃)	Input voltage	110V/220V
Heating temperature	1-5 gears (preheating for 3 minutes, heating temperature 38 to 45 °C)	Output Power	1000W
		Fuse	15A



Function

- ◆ Body Slimming
- ◆ Body Contouring
- ◆ Cellulite Removal
- ◆ Skin Tightening
- ◆ Wrinkle Removal

4 LARGE
HANDLES



Working principle

Infrared RF Vacuum Roller technology combines infrared light, bi-polar radio frequency energy and vacuum, which cause deep heating of the fat cells, their surrounding connective tissue and the underlying dermal collagen fibers. This type of efficient heating and vacuum stimulates the growth of new and better collagen and elastin which results in localized reduction in skin laxity, body volume, and an overall improvement in skin structure and texture.

Advantages

1. Ideal slimming device for circumference reduction.
2. Visible results in as few as 1 treatment session
3. Offers treatment as comfortable as human massag with specially designed vacuum and roller.
4. Delivers optimal results for treating both the deeper tissue and skin's surface.
5. IN 1 technologies combines vacuum+roller+IR laser+cavitation+RF together ,ensure best body shaping and body slimming results.
6. 940nm NearInfrared Laser System

Specification

LCD screen	1) Display screen: 10 inch	RF energy density	Max: 60J/cm
Pulse width	0.5s-7.5s	Laser wavelength	940nm
Negative pressure	1) Absolute value: 80kPa -10kPa (60.8cmHg - 7.6cmHg) 2) Relative value: 20kPa -90kPa (15.2cmHg – 68.4cmHg)	Laser wavelength	940nm
Rev of roller	0-36 rpm	Laser power	MAX 20W
Working mode for roller	4 types	Rated input power	750VA
RF frequency	1MHz/5MHz	Mode of power supply	AC220V-240V~,50Hz/AC100-120V~,60Hz
		Packing method	Wood case

Vacuum+Roller+IR+Cavitation+RF



Function

FACE AND SUBMENTAL



- Skin rejuvenation & skin tightening (Laxity, wrinkles, fine lines)
- Face & neck contouring (submental shaping)
- Skin resurfacing ARMS (BINGO WINGS) & BACK Circumferential reduction
- Skin tightening (Laxity)
- Stretch marks



STOMACH (ABDOMEN & LOVE HANDLES)

- Circumferential reduction
- Body contouring
- Skin tightening (Laxity)
- Stretch marks (red and white)



ARMS (BINGO WINGS) & BACK

- Circumferential reduction
- Skin tightening (Laxity)
- Stretch marks



LEGS, THIGHS & BUTTOCKS

- Circumferential reduction
- Cellulitis treatment
- Body contouring
- Skin tightening

Specification

LCD screen	1) Display screen: 10 inch	RF energy density	Max: 60J/cm
Pulse width	0.5s-7.5s	Laser wavelength	940nm
Negative pressure	1) Absolute value: 80kPa -10kPa (60.8cmHg - 7.6cmHg) 2) Relative value: 20kPa -90kPa (15.2cmHg - 68.4cmHg)	Laser wavelength	940nm
Rev of roller	0-36 rpm	Laser power	MAX 20W
Working mode for roller	4 types	Rated input power	750VA
RF frequency	1MHz/5MHz	Mode of power supply	AC220V-240V~,50Hz/AC100-120V~,60Hz
		Packing method	Wood case

Lift &Tighten (1.5,2.0,3.0,4.5mm) | Tighten &Contouring body(6.0,9.0,13.0mm)

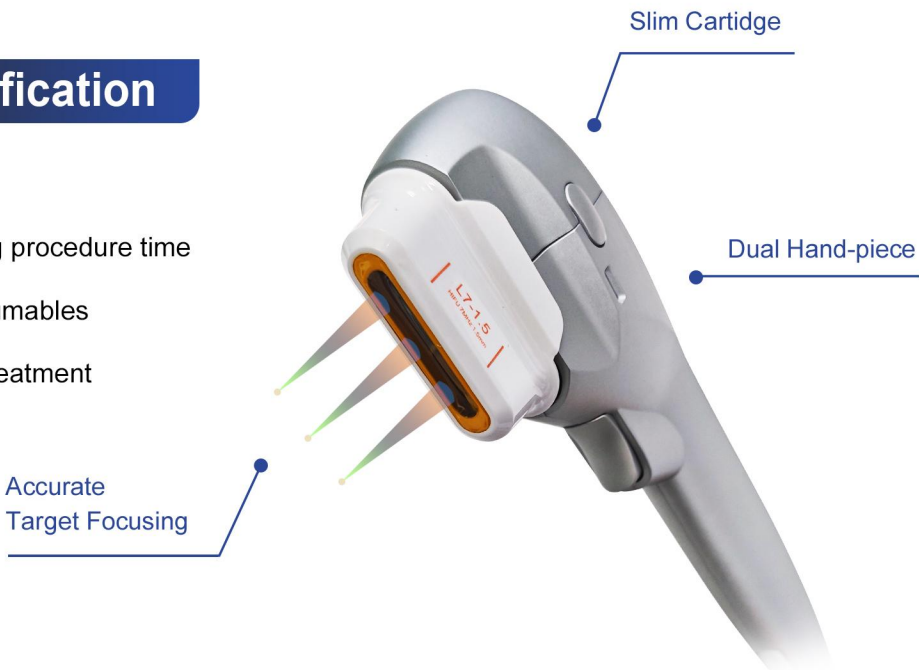


Working principle

The Ultramage HIFU technology with Micro and Macro Focused cartridges, It is a completely non-invasive, non-surgical face lifting and tightening,body tightening and contouring device that provides a customized treatment procedure for the entire face, neck, and body.

Specification

- Shortening procedure time
- Low consumables
- Reliable treatment



Face Treatment

20,000 shots



Body Treatment

30,000 shots



Specification

Target Area	Face, Body, Eyes, Neck/Throat, Lips, Legs/Arms, Vagina	Treatment area	Face,Neck, Body
Operation System	Ultrasonic	Cartridge	7pcs standard
Model Number	TB-7D	Hifu Cartridge option	1.5mm,2.0mm,3.0mm,4.5mm, 6.0mm,9.0mm,13mm
Function	Face lifting Body slimming	Frequency	2Mhz /4 Mhz
Technology	Ultra lift double hifu	Shots	30000 shots/ cartridge
Feature	Skin Tightening, Wrinkle Remover, Anti-Puffiness, Anti-aging, Other	Voltage	110V/220V 50-60Hz
		Dot distance	1-2mm(Set unit 0.1mm)

Taibo HIFU[®]

7D HIFU+RF+VMAX
+VAGINAL HIFU+LIPOSONIC

Eliminate Forehead Wrinkles \ Remove Fine Lines \ Body Slimming
Take Away Nasolabial Folds \ Fat Removal \ Anti Neck Profile



Handle



7D (20000shots)



Vaginal HIFU(10000shots)



VFACIAL MAX(68000shots)



LIPOSONIC(500shots)



RF



Before&After



Before

After

Before

After

Before

After

Can draw any graphics adjustable size and density



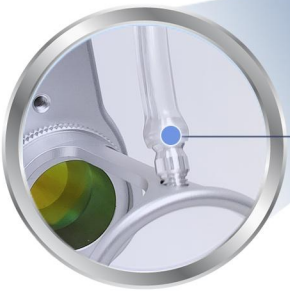
Superiority

- ◆ Scar Removal
- ◆ Wrinkle Removal
- ◆ Skin Rejuvenation
- ◆ Vaginal Tightening
- ◆ Stretch Mark Removal
- ◆ Acne Treatment

360°
ROTATABLE
STAND

360 rotation in three places

High precision scanner
Scan graphics are more accurate and faster



Air outlet
Improve treatment experience and safety



Working principle

Vaginal tightening laser is the use of water uptake 10600nm gold standard laser, using a three-dimensional lattice technology and 360-degree circular emission perfect combination of technology, to produce 50-70 degrees Celsius in a controlled depth of the vaginal mucosa, stimulate lamina propria and myometrium newborn fibroblasts, and damaged collagen fibers, elastic fibers, such as restructuring, so as to achieve vaginal wall thickening and tighteningof the vagina so that firmness, sensitivity enhancement, vagina tight as ever.

Specification

Laser type	Metal tube solid-state laser	Indicating light wavelength	635nm
Laser wavelength	10.6μm±0.1μm	Spot size	≥0.1mm focus
Laser power	1-40W (actual>40W)	Operating Voltage	AC 100-240v~, 50/60Hz
Cooling system	Air cooling	Display size	15.6 inch color touch screen
Light guide system	7 joint light guide arm	Relative working temperature	5℃~40℃
Output mode	Dot matrix mode/pulse mode (long pulse and ultra pulse mode)/private mode/vulvar mode	Relative working humidity	≤80%
		Atmospheric pressure range	700hPa~1060hPa

Multiple modes | Independent operation | Non-interference with each other



Treatment Features

① Non-Ablative SkinTight Applicator



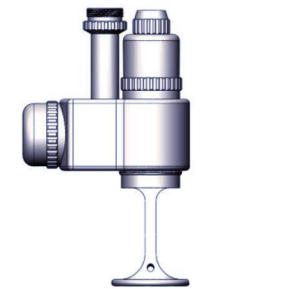
Vagina Applicator Sectinal Vagina Mirror

② Vulva Mode




Cervix Applicator Vulva Applicator

④ Fractional Mode



Deep exfoliation non-exfoliation micro exfoliation

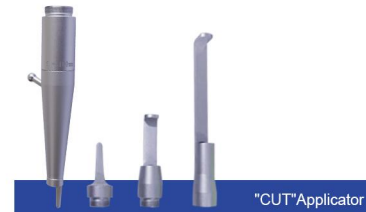③ Super Pulse Mode






Highly intelligent analysis of freckle and acne marks, stretch marks, skin rejuvenation, and scar repair.



Clinical advantages

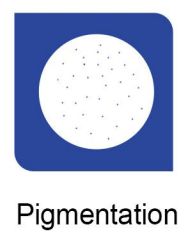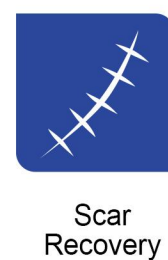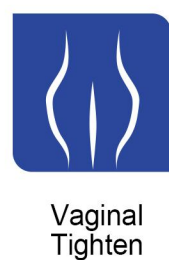- ◆ High energy, ultra short pulse width, less than 2 milliseconds. So the vaporization speed is fast, the thermal damage is minimal, there is no tissue coking, no adhesion formation, and no scar formation. At the same time, the affected area is appropriately heated to stimulate collagen proliferation and rearrangement in the skin, and the clinical healing process is not significantly different or even faster compared to conventional surgical methods.
- ◆ It is possible to accurately grasp the depth of treatment.
- ◆ Less bleeding, as the incision stops bleeding immediately, there is almost no bleeding throughout the entire surgical process, providing a bloodless field of view for precision surgery.
- ◆ The surgery time is short, eliminating the trouble of stopping bleeding, and the wound does not need to be sutured, saving a lot of surgical time.
- ◆ Metal RF tube has uniform beam, stable energy, fine spot, high linearity, and minimal clinical trauma
- ◆ But achieve micro dissection surgery, non dissection



Can draw any graphics adjustable size and density



Superiority

- ◆ Scar Removal
- ◆ Wrinkle Removal
- ◆ Skin Rejuvenation
- ◆ Vaginal Tightening
- ◆ Stretch Mark Removal
- ◆ Acne Treatment



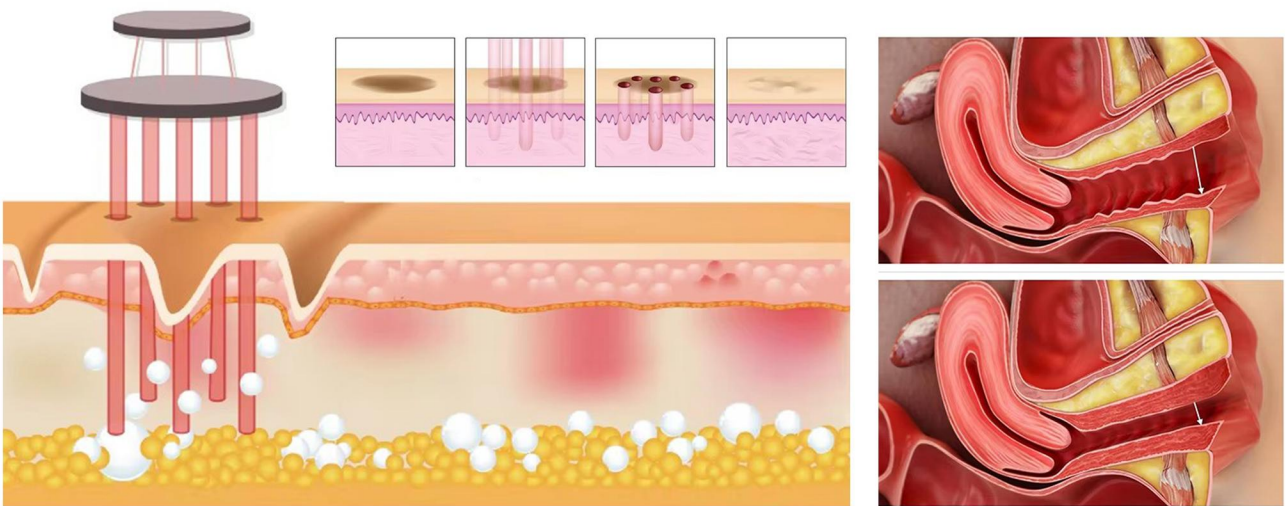
Working principle

Skin rejuvenation >

Fractional laser is a revolutionary progress based on fractional photothermolysis theory and shows the unique advantages in a short time.

Vaginal Treatment principle

Vaginal tightening laser is the use of water uptake 10600nm gold standard laser, using a three-dimensional lattice technology and 360-degree circular emission perfect combination of technology, to produce 50-70 degrees Celsius in a controlled depth of the vaginal mucosa, stimulate lamina propria and myometrium newborn fibroblasts, and damaged collagen fibers, elastic fibers, such as restructuring, so as to achieve vaginal wall thickening and tighteningof the vagina so that firmness, sensitivity enhancement, vagina tight as ever.



Specification

Laser wavelength	10.6μm±0.1μm	Indicating light wavelength	635nm
Laser power	40W	Spot size	≥0.1mm focus
cooling system	Water cooling+heat sink	Operating Voltage	AC 100-120V/AC200-240v~, 50/60Hz
Light guide system	7 joint light guide arm	Display size	10.4 inch color touch screen
Output mode	Dot matrix mode/pulse mode (long pulse and ultra pulse mode)/private mode/vulvar mode	Relative working temperature	5℃~40℃
		Relative working humidity	≤80%
		Atmospheric pressure range	700hPa~1060hPa

Highest Power | Shortest Pulse Width | Versatility for Difficult Cases



7 JOINTS
LASER GUIDE ARM





Tattoos Removal



Freckle Removal



Pigment Removal



Skin Rejuvenation



Hair Removal



Vascular Removal



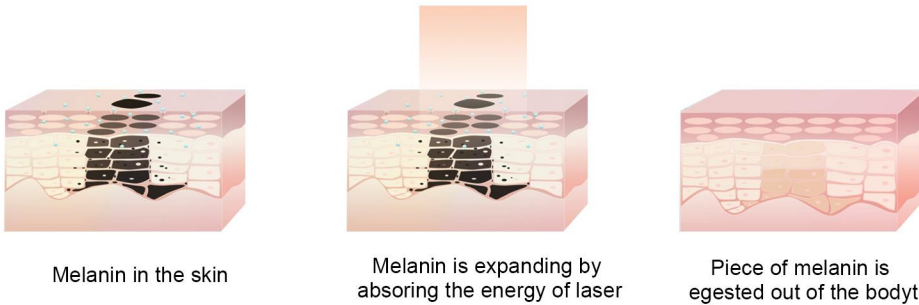
Onychomycosis



Nevus Of Ota Removal

Working principle

The PicoSecond Laser uses ultra-short pulses (one trillionth of a second in length) to hit the melanin with a great pressure, the melanin shatters into tiny dust-like particles. Because the particles are so small, they are more readily absorbed and eliminated by the body. It is a quick and easy non-surgical ,non-invasive laser skin treatmten for the body.



Specification

Output mode	Pico Laser Q-switch Pulsed output	Work mode	Continually work with intervals
Wavelength	1064nm, 532nm, 755nm	Spot Size	2-10mm
Pulse Duration	6ns	Output energy	1064nm :100-2000mj
Frequency	1-10Hz adjustable		755nm :100-2000mj
Single Pulse energy	100mj-2000mj		532nm :50-1000mj
Safety distance	NOHD 1km Operator must wear goggles	Wastage power	1000VA Fuse 220V, 10A
Cooling system	Water cooling+air heat exchange		

Highest Power | Shortest Pulse Width | Versatility for Difficult Cases



Products advantages



Professional & Effective
Freckle Tattoo Removal



Adjustable Energy
& Spot size



Full Solution For
Variety Skin Problems



Laser guide arm

10
Sizes



FAST AND POWERFUL



TIME IS SHORT
LESS DAMAGE



ABSORB WELL
RECOVER QUICKLY

7 joints articulated

Steelyard Weight

Automatically raise light guide arm


Specification

Output mode	Pico Laser Q-switch Pulsed output	Work mode	Continually work with intervals
Wavelength	1064nm, 532nm, 755nm	Spot Size	2-10mm
Pulse Duration	6ns	Output energy	1064nm :100-2000mj
Frequency	1-10Hz adjustable		755nm :100-2000mj
Single Pulse energy	100mj-2000mj		532nm :50-1000mj
Safety distance	NOHD 1km Operator must wear goggles	Peak Power	3000W
Cooling system	Water cooling+air heat exchange	Voltage	220V or 110v ,50hz
Weight	105kg	Dimensions	106x53x106cm

Highest Power | Shortest Pulse Width | Versatility for Difficult Cases



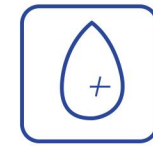Versatility for Difficult Cases



Tattoos Removal



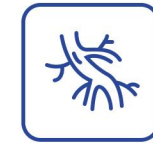
Freckle Removal



Pigment Removal



Hair Removal




Vascular Removal



Onychomycosis



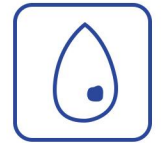Skin Rejuvenation



Nevus of OTA Removal

Detail Display



Break-through 5 working modes

PICOTECH

Pigments, Tattoos, Birthmarks, Skin Rejuvenation

PTP

Melasma, Moles, Tattoo, Hyperpigmentation.

MULTI PULSE

Chloasma, Stubborn Pigments.

LONG PULSE

Vascular, Scars, Acne, Rosacea, Diffuse Redness, Onychomycosis, Rejuvenation, Firming, Whitening, Shrinking Pores, Remove Fine Wrinkles, Hair Removal, Improve Veins.

6+1 980+1470NM DIODE LASER MACHINE+OPTICAL FIBER

Non-invasive comfortable | 6+1 Functions



Working principle

Vascular removal

The instrument uses a 980nm wavelength laser to transmit laser energy to the blood vessel part through an optical fiber. According to the principle of selective photothermal, hemoglobin absorbs laser energy and converts it into small molecular clusters and is absorbed by tissues.

Pain physiotherapy

Focal luminescence produces thermal stimulation, uses the biological effect of laser to act on the human body, enhances capillary permeability, and increases ATP production , activates healthier cells or tissues, achieves analgesia, speeds up tissue repair and healing.

Onychomycosis removal

Laser treatment is a new type of treatment that uses the principle of selective light and heat to kill fungi without damaging normal tissues and has no side effects. With anti-fungal ointment, the effect is better.

Skin rejuvenation &anti-inflammatory

Due to the effect of improving the phagocytic function of white blood cells, it can affect the activity of enzymes, regulate the immune function of the body, and achieve anti-inflammatory The purpose of reducing swelling and speeding up the tissue repair process.

Eczema and herpes

The laser beam generated by the semiconductor laser continuously and directly irradiates the patient's skin lesions, and its laser energy can be absorbed by the tissue and converted into bioenergy, induce or activate macrophages and lymphocytes, improve the effects of specific and non-specific immunity.

Endolifting,Lipolysis surgery,EVLT surgery

The instrument mainly acts on deep fat, superficial fat, and transfers the energy directly to the fat cells for uniform heating. During the heating process, the connective tissue and fat cell structure can be changed by controlling the heat, and the adipose tissue has a photo thermal effect (so that the fat is dissolved). Meanwhile, the photodynamic effect (separating the fat cells from the normal tissue) decomposes the fat cells to make them evenly liquefied, and the fat liquid is excreted through the ultra-fine positioning needle, which fundamentally reduces the number of fat cells, effectively avoids postoperative rebound.

Specification

Output	Fiber-optic coupling	Frequency	1-30Hz
Aiming	Beam 650nm	Operation interface	8.0 inch
Laser wavelength	980nm	Mode	Pulse mode, Continuous mode
Net/Gross weight	5KG/10KG	Cooling	Air cooling
Output Power	1-60W(options)	Operate mode	Touch screen Machine
Input AC:	100-240V、 50/60Hz	Dimensions	360mm*340mm*220mm
Pulse width	15ms—100ms	language	English or OEM
Length of Fiber	2m,3m	Flight case Dimensions	460mm*440mm*340mm



M8 MICRONEEDLE RF MACHINE

Reduce wrinkles | Skin Tightening | Skin Whitening | Pore Reduction



- Reduce wrinkles
- Remove acne & scars
- Remove stretch marks
- Shrink pores
- Lift and tighten skin
- Acne removal



Treatment handle

7mm



- 12p: For small or smaller treatment areas.
- 24p: For face and neck.
- 40p: For use on body parts for larger and deepertissue treatments.
- Nano crystal head:The introduction enhancesand tightens.

Working principle

Fractional RF Technology allows effectively treatment of deeper dermis, epidermis,and upper dermis. RF energy is emitted at the end of fractional RF microneedle; this will treat target area - deeper dermis - without damage of epidermis. The fractional RF is useful for minimal invasive epidermal rejuvenation. RF energy does not target pigment-making cells in the skin. Cures and improves wrinkle, scar, skin re-juvenation,and skin tightening safely and efficiently with minimized PIH, using RF energy.

Specification

Energy	face	0-200W
	body	24/40 needles: 0-300w
		12 needles/NANO: 0-200W
Needle depth	face	0.5-2.5mm
	body	24/40 needles:0.5-7mm
		12 needles/NANO: 0.5-2.5mm
Interface size	12.1-inch	
Input voltage AC	110V~220V~50Hz/60Hz	
Treatment probe	12p 24p 40p nanocrystal head	

4 IN 1 HYDRO DERMABRASION MACHINE

Skin Cleaning | Skin Revitalizer | Exfoliators | Blackhead Removal



SUPER SUCTION

70kpa

Super suction, take away stubborn blackheads

The suction power can reach up to 100kpa, which can strongly absorb dirt and blackheads in the skin and eliminate dead skin.

ordinary suction

others

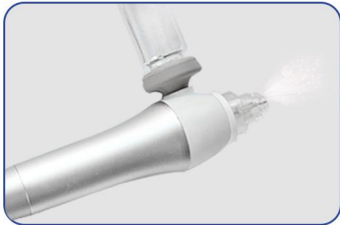
Super suction

our machine






Working principle




ION SPRAY GUN

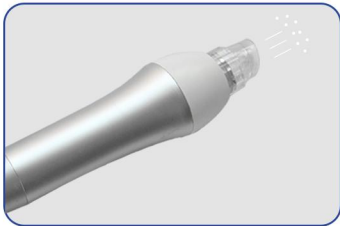The ion generator continuously decomposes hydrogen and oxygen ions from the water, and the spray gun atomizes the original liquid into nutrient water molecules. With the hydrogen chloride ion medium, the skin absorbs nutrients more easily and the skin is more resistant to oxidation and shines.



BUBBLE PEN

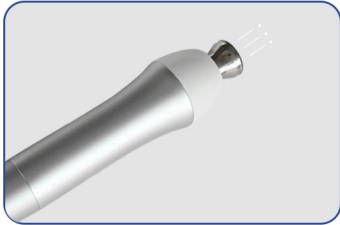Improve stain, color sink, shrink pores, dark skin, yellowish, blackheads, deep cleansing, toning, enhance skin elasticity luster, tender skin, etc



DIAMOND HEAD

removes skin debris, scars, blemishes, wrinkles and uneven pigmentation of the skin, help to replenish natural nutrients and restore cell activity promoting healthy skin.



RF HEAD

Improves fine lines, face lift, heck skin lifting, skin tightening,smooth wrinkles.

Specification

Vacuum Power	70Kpa Max	Operation Control	Control Panel
Technology	Water dermabrasion & Diamond dermabrasion	Electrical	AC220V/50Hz AC110V/60Hz
Technology	1) Water facial / Waterdermabrasion	ion spray	10ml
	2) Microdermabrasion	RF power	10w
Handpieces	1)Waterdermabrasion handle with 8 tips	Machine size	300mm*380mm*360mm
	2) Microdermabrasion handle with 9 tips	Machine weight	10kg(Not counting packaging)
	3) Ion spray gun	Package size	460mm*460mm*460mm
	4) RF handle		